


# LumaCare

Illumination for Photodynamic Therapy  
plus PhotoTherapeutic and  
Cosmetic Applications

A simple non-coherent  
light source delivering all  
frequencies from one device.

FDA Approved  
CE MDD Approved

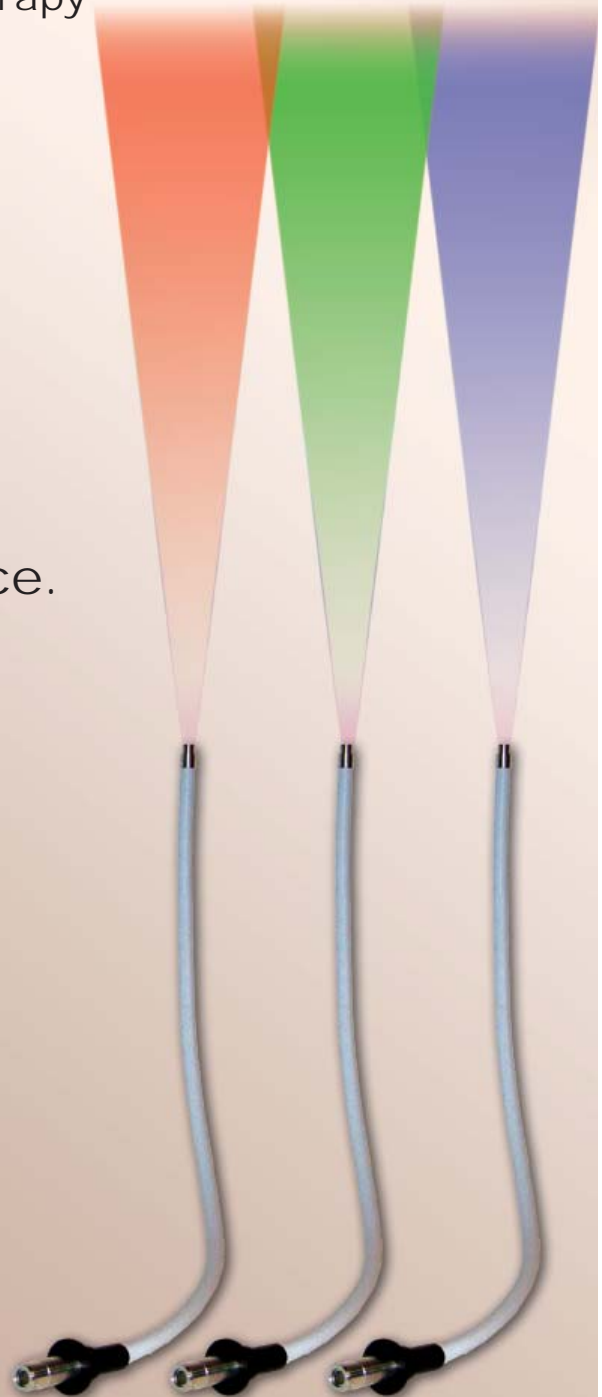

# Single Solution

LumaCare® has developed an illumination source for Photodynamic Therapy that is simple, elegant, safe and cost effective. Our Patented technologies utilize a pure white, full spectrum light source to create the full spectrum of potential light frequencies. Then we filter and focus the light source with a fiber optic probe to deliver protocol specific frequencies and beam patterns required for most PDT protocols.

This patented technology has two key components; the Light Source produces the full spectrum of visible light, and the Fiber Optic Probe(s) filter and focus the light to specific frequencies, beam patterns and power levels.

## Probes

LumaCare LC-122 offers interchangeable, protocol specific probes that are easily connected with a simple interlocking connection.

Doctors and Staff can now implement PDT in their existing practice in a cost effective and safe manner. The cost effective Fiber Optic Probe(s) are offered in a wide range of frequencies or power requirements and are easily changed in seconds. This allows Doctors to utilize many PDT drugs, compounds and protocols without purchasing numerous lasers or drug specific single use light sources.

## The Light Source

The Light Source is simple, consisting of the non-coherent light source capable of producing almost the entire spectrum of visible light. Having the entire visible light spectrum available allows the LC-122 to produce almost any frequency of light for a wide variety of PDT protocols. The LC-122 can be taken into multiple treatment rooms to increase the number of patients treated. The LumaCare light source is not a laser or LED array. LumaCare eliminates the high cost of ownership and removes the LED array / laser single frequency impediment.

	Applications*
■ <b>Medical Treatments*:</b>	<b>Acne, Actinic Keratoses (AK), Basal Cell Carcinoma (BCC), Psoriasis, Wound Healing, Dermatological Lesions, Wrinkle Reduction and Cosmetic Protocols</b>
■ <b>Direct Replacement for:</b>	<b>Aktelite™ (all models), Clearlight™ (all models), Diomed™ (all PDT lasers), DUSA Blu-U™, Lightwave Delux and Omnilux™</b>
■ <b>Approved* Protocols:</b>	<b>ALA, 5-ALA, Levulan™, Metvix™, Photofrin™, Visudune™</b>

\* Where Authorized

## Specifications

Voltage:	85VAC – 260VAC	Lamp type:	Specially Designed Quartz Halogen (mixed) lamp
Current:	1.6A -3.2A	Operating Temperature:	+10° to +40°C
Frequency:	47Hz – 63Hz	Storage Temperature:	-30° to +60°C
Cooling:	Forced Air Cooling ~50m³/hr	Relative Humidity:	5% to 95%

LumaCare is a patented medical device covered by one ore more of the following International and USPTO Patents: 5,849,027, 6,187,030, EP 1124613 B1.



**LumaCare Medical Group**  
USA Phone: +1 (949) 644-0126

[www.LumaCare.com](http://www.LumaCare.com)

A Division of MBG Technologies,  
[info@LumaCare.com](mailto:info@LumaCare.com)